EndoTherapy

EndoTherapy

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- _ Wiederaufbereitung
- _ Reprocessing
- _ Reconditionnement
- _ Reacondicionamiento
- _ Ricondizionamento
- _ Opnieuw gebruiksklaar maken
- _ Reprocessamento
- _ Upparbetning
- _ Rensning
- _ Επανεξεργασία
- _ Yeniden hazırlama

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Recommendations for the reprocessing of endoscopic instruments

In order to guarantee the functional reliability and a long service life, reusable endoscopic instruments must be conditioned after each examination in accordance with the applicable guidelines and standards. These recommendations are based on the E.S.G.E.N.A "Protocol for reprocessing endoscopy accessories", "Hygiene requirements for the conditioning of flexible endoscopes and endoscopic additional instruments" of the Robert-Koch Institut (RKI) and on DIN EN ISO 17664.

Please note that the operator is solely responsible for carrying out and monitoring the correct and proper reprocessing.

Products

The recommendations published here contain important information for the reliable and effective reprocessing of medwork® products and apply in conjunction with the respective instructions for use of the products.

All reusable medwork® products which are marked with this symbol "134 °C" can be supplied for reprocessing. Products which are intended to only be used once (☉-symbol) must never be reprocessed to the reprocessing process.

Instruments which penetrate tissue or hollow organs must be sterile. Instruments which contact intact mucous membrane must be disinfected. In the case of instruments with joints (biopsy forceps, foreign body forceps), it must be ensured that the branches are open during the complete cleaning process.

It must particularly be ensured that standardised and validated processes are applied. If reusable medwork® products are used for patients suffering from a variant of the Creutzfeldt-Jakob disease, these must not be conditioned after use. The particular recommendations of the RKI are applicable here.

The effectiveness of the cleaning and disinfectant agent gigasept AF forte® (Schülke & Mayr GmbH) has been confirmed in extensive tests and examinations for the conditioning of medwork® products. In the case of using other cleaning and disinfectant agents, always ensure that these are suitable for this use.

Premises

Only perform the reprocessing in suitable rooms which are specially equipped for this purpose and which comply with the requirements of the Committee for Hospital Hygiene and Infection Prevention at the Robert-Koch-Institut (RKI).

Personnel

Please note that endoscopic additional instruments must only be conditioned by trained persons (Category IV) who have the necessary training, level of knowledge and experience for this. You should wear face protection, an operating gown impervious to liquids and cut-resistant gloves for protection against possible contamination with pathogenic germs and against potentially hazardous cleaning chemicals.

Step 1: Cleaning

Precleaning:

We recommend here:

Equipment: cleaning agent, brush, disposable cloth or sponge, mains water with at least drinking water quality

Cleaning agent:

- Cleaning and disinfection agent
gigasept AF forte® (Schülke & Mayr GmbH; application concentration 2%; application time 15 minutes)
- Surfactant-based cleaning agent
Bodedex forte® (Bode Chemie Hamburg; application concentration 0.5% - 1%; application time 5-10 minutes depending on degree of soiling)
- Enzymatic cleaning agent
gigazyme® (Schülke & Mayr GmbH; application concentration 0.5% - 1%; application time 5-10 minutes depending on degree of soiling)

Clean the used instrument as quickly as possible after the endoscopy with a disposable cloth or sponge (e.g. MIS1 Cleanscope) in order to minimise as much as possible drying on of organic material and chemical residues and contamination of the environment. For the protection of personnel and the environment, only transport contaminated instruments in closed containers from the examination room to the reprocessing room.

If possible place the instruments immediately after use in a surface-active, non-fixing cleaning solution with strict observance of the concentrations, interactions and application times stated by the manufacturer. Only use cleaning solutions for this whose suitability and effectiveness have been documented and declared accordingly in the expert reports of the preparation and equipment manufacturers. Do not use any solutions containing aldehydes as these can fix proteins and thus adversely affect the success of the cleaning. We discourage the use of alkaline cleaning agents as these adversely affect the material quality of the products and thus shorten the product service life. Use water of at least drinking water

quality for creating the solution. First remove heavy external soiling with a soft cloth, sponge or soft brush soaked in cleaning solution. The cleaning solution must be replaced at least daily and immediately in the case of visible contamination.

Dismantle the instrument as far as possible and carry out all other cleaning steps underneath the liquid level in order to prevent spraying of contaminated liquid. Rinse afterwards with a syringe in each case of at least 20 ml cleaning solution through all lumen and hollow spaces in order to remove secretions and tissue residues. In doing so, ensure that all lumen and hollow spaces are completely flushed and no air inclusions form.

ATTENTION! Special Instructions

medwork® products with rigid handle must not be dismantled. Rinse the lumen sufficiently using the luer lock attachment on the handle.

Ultrasonic Cleaning

We recommend here

Equipment: ultrasonic cleaner (operating frequency 30-50 Hz); cleaning agent; mains water of at least drinking water quality

- Cleaning and disinfection agent
gigasept AF forte® (Schülke & Mayr GmbH; application concentration 2%; application time 5 minutes)
- Surfactant-based cleaning agent
Bodedex forte® (Bode Chemie Hamburg; application concentration 0.5% - 1%; application time 5-10 minutes)

Take the instruments out of the cleaning solution and load the basket of the ultrasonic bath with them. Ensure that its dimensions allow complete submersion of the components in the cleaning solution. Also note here that too heavily loaded baskets can result in dead areas which are not reached by the ultrasonic waves.

ATTENTION! Special Instructions

If you use enzymatic cleaning solutions for the ultrasonic cleaning, ensure that the cleaning agents are verifiably suitable for use in the ultrasonic bath and are not foaming. The temperature range recommended by the manufacturer must also be strictly complied with. Ensure that the stated maximum temperature is not exceeded. The temperature in the ultrasonic bath should be checked and regulated by the machine for this. It is recommended to produce a cleaning report. The cleaning and disinfection solution in the ultrasonic bath must be replaced at least daily and immediately in the case of visible contamination.

- _ Insert the instrument with a diameter of at least 20 cm into the basket.
- _ Fill all lumen and hollow spaces with the cleaning solution and avoid bubbles.
- _ Flush the lumen of the instrument after the operation in the ultrasonic bath with the solution until all soiling has been completely eliminated.
- _ Remove the instrument from the ultrasonic cleaner and repeat the flushing process with water of drinking water quality until all soiling and cleaning agent residues are completely removed.

Step 2: Disinfection

We recommend here

Equipment: disinfectant, mains water with drinking water quality

- Cleaning and disinfection agent
gigasept AF forte® (Schülke & Mayr GmbH; application concentration 2%; application time 15 minutes)
- Aldehyde-based disinfectant
forte® extra; (Bode Chemie Hamburg; application concentration 4%; application time 60 minutes)

Place the cleaned instruments in a pan with disinfectant solution.

- _ Afterwards, fill all lumen free of air bubbles with disinfectant solution.
- _ Seal the pan with the associated lid.
- _ Strictly observe the concentration and application time of the disinfectant in accordance with the manufacturer information.
- _ Remove the instruments from the disinfectant solution using new disposable gloves.

Step 3: Rinsing / Washing machine treatment

Rinsing

- _ Place the cleaned instrument in demineralised water of drinking water quality. Replace the water after each cleaning cycle.
- _ Rinse the instrument with water of drinking water quality and afterwards rinse the inner lumen several times with a syringe sufficiently long to remove all residues.
- _ Remove the instrument from the water.

Washing machine treatment

- _ Load the thoroughly cleaned instruments in the washing machine containers provided for this. Observe the manufacturer information.
- _ Place handles, wires and rolled up catheters in the special basket provided for this and connect the instrument with lumen to the provided connections in accordance with the manufacturer information in order to ensure complete and thorough rinsing of all hollow spaces.
- _ Do not close the washing machine until after you have removed your gloves.
- _ Select the suitable washing cycle and start the program.
- _ After completion of the washing program, ensure that all cycle steps have been performed and that the parameters have been complied with in accordance with the manufacturer information.
- _ Open the washing machine and remove the instrument.

Step 4: Drying

- _ The surface of the catheter can be dried with a fluff-free cloth. Each lumen must be completely dried using medical compressed air.
- _ Suspend instruments with spiral casing (e.g. biopsy forceps) freely in order to support the drying process.

Check the instruments for function, damage and wear afterwards. Pay particular attention here to functional movement of the components and check the completeness of all components for instruments which can be dismantled.

Attention! Do not use any products which show kinks, cracks, fractures, discoloration, surface changes or similar.

Step 5: Sterilisation

- _ Assemble dismantled instruments completely.
- _ In order to improve the sliding ability and functionality of the components, we recommend suitable care agents (e.g. No. 460044 – Medicon/Tuttlingen).
- _ Pack the instruments in suitable, validated sterile packaging.
- _ Select the suitable sterilisation method for the instrument (e.g. pre-vacuum, 134 °C, holding time: 5 minutes or corresponding cycles) in accordance with the manufacturer information and national guidelines. All reusable medwork® products except the oesophagus bougies (DIL2-D2-...) can be sterilised in an autoclave.
- _ We do not recommend coiling more tightly than 20 cm in order to prevent damage to the instrument during the sterilisation, transport and storage.
- _ After completion of the sterilisation process, ensure that all cycle steps have been performed and that the parameters have been complied with in accordance with the manufacturer information.
- _ Check the sterile packaging for damage and sterilisation indicators.

Step 6: Storage

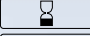
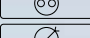
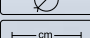
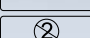
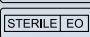
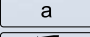

- _ All medwork® products must be stored in a dry place, protected from light at +15 °C to +25 °C and at relative humidity of 35% to 50%.
- _ Keep all usage instructions easily accessible.
- _ Before use, check the sterilisation date on the packaging as sterile products must be used by a certain date.
- _ Do NOT use the product if the sterile packaging shows cracks or perforations, the seal is damaged or moisture has penetrated.

Further information about reprocessing can be found at:

- _ www.esgena.org
- _ www.rki.de

Recommended guide parameters for the autoclave cycle

Loading	Standard
Process	Pre-vacuum
Temperature	134 °C
Sterilisation pressure	2.5 bar
Sterilisation time	5 minutes

	Manufacturing date
	Use Before
	Attention, see instructions for use
	Lot number
	Number of lumens
	Exterior diameter
	Usable Length
	For single use only
	Catalogue number
	Ethylene oxide sterilized
	Can be sterilized up to 134 °C
	Store in a dry, dark place between 15°C and 25°C
	Variable product data
	Max. guiding wire diameter
	Non-sterile product

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