

medw|ork®




icus



The **new design** of **LCUS** injection needles reflects the increasing demand for ergonomics and functionality in clinical everyday life. Added to this they provide what really matters: reliability and safety.

Safe needle guidance

The shielded needle prevents damage to the working channel. Even if the endoscope is in an angled position, the needle always remains centred in the catheter.

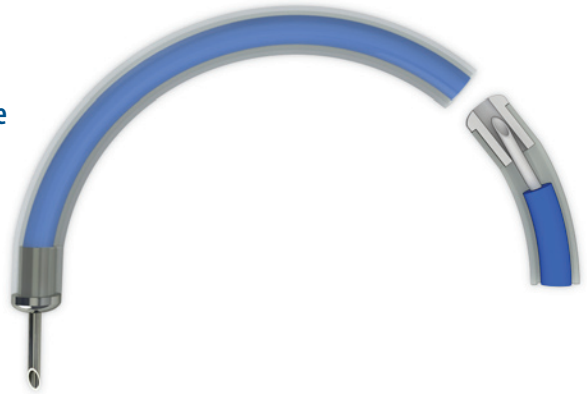

Easy puncture

The proven lancet tip makes it easier for you to puncture the tissue smoothly.



Tactile feedback and ergonomic design

Feel the distinct engagement of the needle in the retracted and extended positions.



Disposable injection needles

Product Code	ø Needle	Needle-lift	øTube	Length	Unit
INJ1-A1-05-5-18-180	0,5 mm	5 mm	1,8 mm	180 cm	10
INJ1-A1-07-4-23-230	0,7 mm	4 mm	2,3 mm	230 cm	10
INJ1-A1-07-5-18-180	0,7 mm	5 mm	1,8 mm	180 cm	10
INJ1-A1-07-5-23-180	0,7 mm	5 mm	2,3 mm	180 cm	10
INJ1-A1-07-5-23-230	0,7 mm	5 mm	2,3 mm	230 cm	10
INJ1-A1-07-8-23-230	0,7 mm	8 mm	2,3 mm	230 cm	10
INJ1-A1-09-5-23-180	0,9 mm	5 mm	2,3 mm	180 cm	10
INJ1-A1-09-5-23-230	0,9 mm	5 mm	2,3 mm	230 cm	10

Disposable injection needles for balloon enteroscopy

Product Code	ø Needle	Needle-lift	øTube	Length	Unit
INJ1-A1-07-5-18-260	0,7 mm	5 mm	1,8 mm	260 cm	10
INJ1-A1-07-5-23-260	0,7 mm	5 mm	2,3 mm	260 cm	10

